

TREATMENT - UNIT DOSE MEDICATION - SIMPLE TROLLEY 100 SERIES 1817 W - 10 PATIENTS

DIMENSIONS 640 x 480 mm

ADEQUATE HEIGHT, ERGONOMICS 1005 mm (work surface – floor).

WORK SURFACE

- Polyethylene, PE. Asepsis, hygiene, easy-to-clean.
- Withstands alcohol, detergents, disinfectants, chlorhexidine.
- Protection on all 4 sides. Pleasant to the touch, rounded shapes.

MONOBLOC STRUCTURE

- Stainless steel AISI-304. Assembly with rivets and TIG welding.

ISO DRAWER FRAMEWORK



- Total opening.
- Steel ball-bearing slides with self-closing system.
- Front stop, it protects the slide against strong impacts when closing.
- Load capacity per drawer: 20 kg.
- Interior frame ISO 600 x 400 mm.
- Polyethylene drawer front part, thickness 10 mm, 6 colours (Ref. G, Y, W, B, F, R).
- Ergonomic handle.
- Card holder with white cards, removable.

ISO DRAWERS 600 x 400 mm

- Heights 50, 100, 200 mm. Interchangeable, ABS.
- ISO norm 3394.

INTERCHANGEABLE CASSETTES, 5 PATIENTS

- Compatible with pharmacy ISO furniture
- Polypropylene and aluminium, grey colour
- Cassette lock, epoxy painted steel
- Locked cassette, card upper part hidden, lower part visible.
- Cassette total weight 4,6 kg.

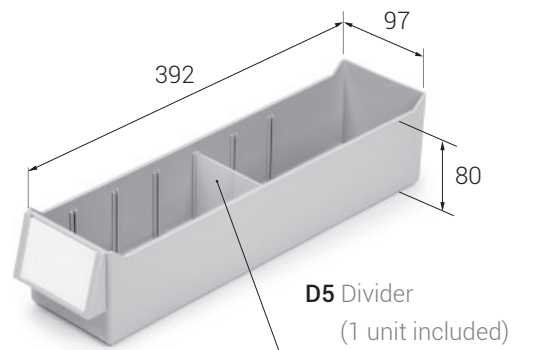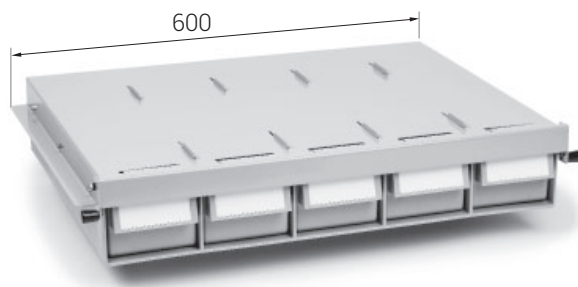

- Capacity: 3 litres.
- ISO norm 7810 identification card, international standard
- Appropriate for intravenous therapy medication, re-packaged medication in unit doses, pills, prefilled syringes, syrups, etc.

KEY LOCKING

- Total locking of ISO drawers and cassettes for medication.
- Prevents from opening ISO drawers, the medication bins and the cassettes extraction.

EXTENSIVE MOBILITY TROLLEYS

- Elevated push handle. Ergonomics, height 1045 mm.
- Double wheels Ø 125 mm, 2 antistatic, 2 with brake.



1817 W - 10 PATIENTS

108M6 W.26(RA4.A5).27RG.(2)32L.33R.50

ACCESSORIES

26(RA4.A5) BINS DISPENSERS

27RG LAPTOP ADJUSTABLE BASE

(2)32L SIDE SUPPORTS

33R DISPOSABLE CONTAINER SUPPORT

50 13 LITRES BIN